

Medical Solutions' competences

Disposable Breathing Bags



One of the largest breathing bag producers in the world


- Widely accepted product design
- Automated production system
- Consistent quality
- Latex-free material
- Conforms to ISO 5362
- 1/2 L Size
- 1 L Size
- 2L Size
- 3L Size
- Available with rigid or soft ISO connector
- Available without connector

Nolato offers a wide range of high quality breathing bags manufactured in a cost efficient medical device environment.

Nolato is one of the largest producers in the world for anesthesia breathing bags and selected as a partner by the market leaders in the anesthesia machinery and breathing circuit markets.

Nolato's breathing bags are made from synthetic latex material. We offer a complete range of latex-free bags in sizes ranging from 0,5 ltr, 1,0 ltr, 2,0 ltr and 3,0 ltr with different types of connectors in compliance with ISO regulations. Breathing bags are available either as disposable or reusable.

Nolato has a modern and highly automated production line for breathing bags, which ensures high and consistent quality as well as high capacity.

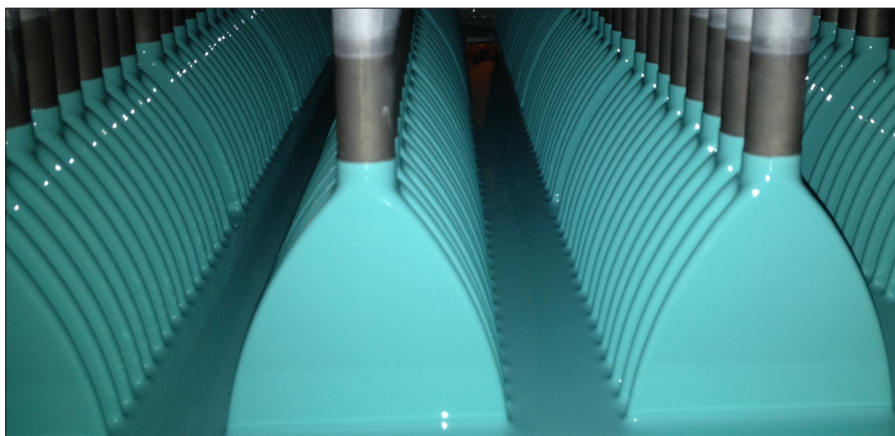



1/2 L



1 L



2 L




3 L



Soft TPE connector



Rigid PP connector

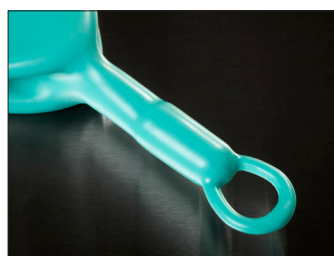

Loop tail



Extended tail 1/2 L bags

For more information, RFQ, orders etc. please contact:

Carolina Kjellman
Product Manager
Phone + 46 729 83 91 60
E-mail: carolina.kjellman@nolato.com

Breathing bags are manufactured in accordance with GMP, ISO 13485, ISO 9001 and ISO 14001 requirements. Each individual breathing bag is quality inspected by Nolato.

Product standard fulfillment:
ISO 5362 (Anaesthetic reservoir bags)
ISO 5356-1 (Conical connectors)

Biocompatibility and raw materials:

Nolato's range of breathing bags are made from a latex-free polychloroprene rubber compound biocompatibility tested in accordance with ISO 10993:

- Skin sensitization test
- Intracutaneous reactivity
- Cytotoxicity by elution



Nolato MediTor, SE-269 78 Torekov, Sweden | Phone +46 431 44 22 60 | www.nolato.com

Nolato MediTor factsheet#7, Rev 3 February 2025. Nolato assumes no responsibility or liability for any errors in the product information.